

PURSUIT[®] XL

ENGINEERED FOR POWERFUL
PERFORMANCE AND
EXCEPTIONAL STYLE

PURSUIT[®] XL

The **Pursuit[®] XL** is an exciting blend of power and precision that features full suspension, low-profile tires and a hydraulic sealed brake system. Standard luxury touches like a wrap-around delta tiller and LED lighting make it the ultimate outdoor scooter.



Delta tiller for optimal handling



LED lighting package



Large 13" pneumatic tires



Full suspension



Pursuit® XL

Model number:	SC714
Weight capacity:	400 lbs.
Maximum speed ¹ :	Up to 8.25 mph
Ground clearance ² :	5.25"
Obstacle climbing ability:	4.72"
Turning radius:	82.5"
Overall length ² :	56.25"
Overall width ² :	27.25"
Front tires:	13" pneumatic
Rear tires:	13" pneumatic
Seat to deck:	21"-22.5"
Suspension:	Front & Rear
Range per charge ¹ :	31.83 miles with Group 31 (110 AH) 28.94 miles with Group 27 (100 AH) 21.70 miles with Group 24 (70 AH)
Total weight without batteries:	205.5 lbs. (main frame)
Standard seating:	Type: CRS reclining high-back seat with headrest Weight: 51 lbs. Material: Black Vinyl (standard) Dimensions: width - 18" (overall width) depth - 18" (effective depth)
Drive system:	Rear-wheel drive, sealed transaxle, 24-volt DC motor
Dual braking system:	Regenerative and electromechanical
Battery requirements ³ :	Two 12-volt, deep cycle Size: Group 31 (110 AH) Group 27 (100 AH) Group 24 (70 AH) Weight ⁴ : Group 24: 53.5 lbs. each; Group 27: 63.5 lbs. Group 31: 70 lbs.
Battery charger:	8-amp, off-board

Colors:



Candy Apple Red



Viper Blue



Silver

¹Varies with user weight, terrain type, battery amp-hour (AH), battery charge, battery condition and tyre condition.

²Due to manufacturing tolerances and continual product improvement, this specification can be subject to a variance of (+ or - 3%).

³AGM or Gel-Cell type required.

⁴Battery weight may vary based on manufacturer.

NOTE: All specifications subject to change without notice.

Features


- Infinitely-adjustable tiller angle with conveniently located handle
- Powerful drivetrain, features 180 amp microprocessor based controller
- CRS high-back seat w/ headrest and sliders
- Full LED lighting package that includes directional signals, automatic brake lighting and built-in, high-intensity LED front light
- Delta tiller with wraparound handles and quick-adjust lever
- Large 13" low profile pneumatic tires for excellent outdoor performance
- Front and rear suspension
- Easily accessible tie down points (for transport of unoccupied scooter)
- Front basket standard
- Rearview mirror standard
- Hydraulic front disc brakes

Options

- Rear basket
- Safety flag
- Durable weather cover
- Special measurement seats available

Pride™ Lifts & Ramps

- Single Fold Ramps
- Solid Ramps

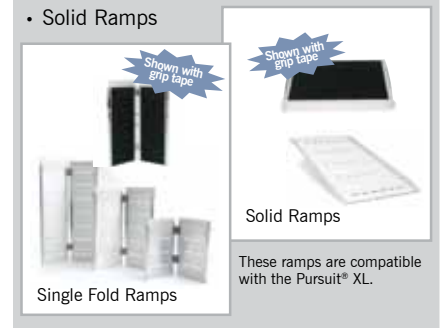

The information contained herein is correct at the time of publication; we reserve the right to alter specifications without prior notice. Speed & range vary with user weight, terrain type, battery charge, battery condition and tire pressure. Due to manufacturing tolerances and continual product improvement, length, width, turning radius and ground clearance can be subject to variance of + or - 3%.

800-800-8586 (US) • 888-570-1113 (Canada) • www.pridemobility.com

